


Space0AR™ Hydrogel

Prostate-rectum spacing



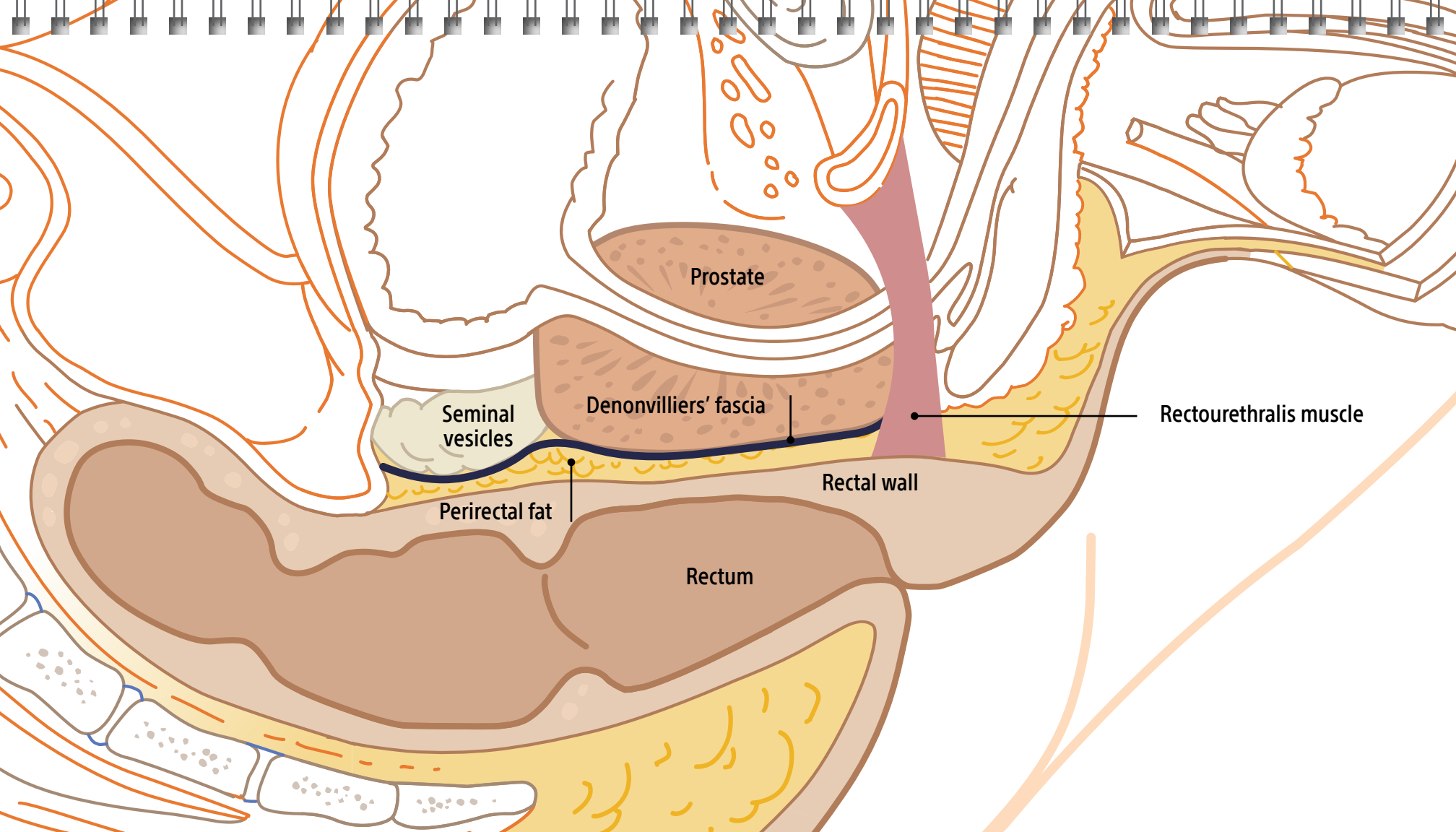

Prostate

Seminal
vesicles


Denonvilliers' fascia

Rectourethralis muscle

Rectal wall

Perirectal fat


Rectum

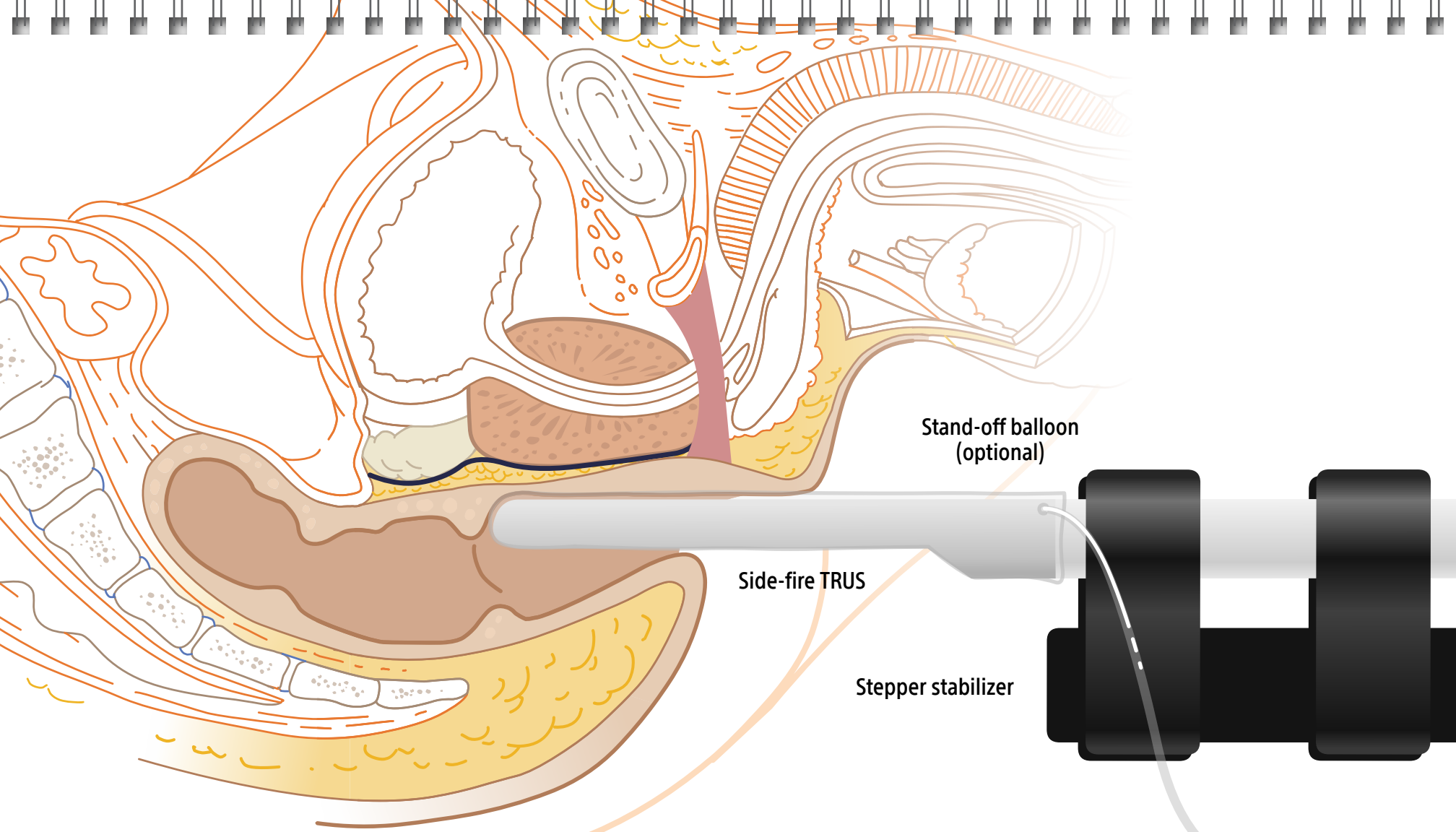

Stand-off balloon
(optional)

Side-fire TRUS

Stepper stabilizer

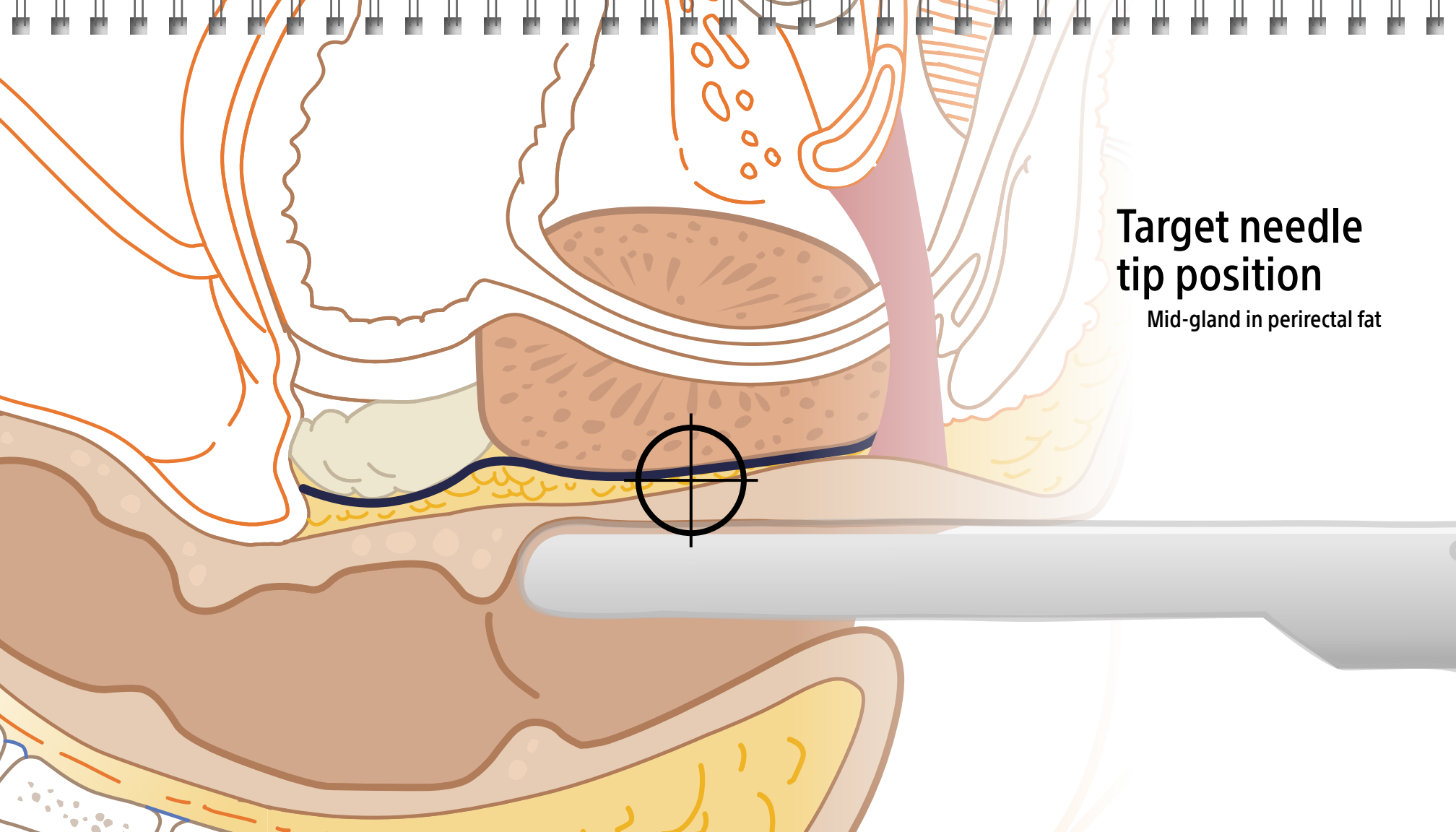


**Target needle
tip position**


Mid-gland in perirectal fat

Ultrasound - Sagittal view

**Target needle
tip position**




Mid-gland in perirectal fat

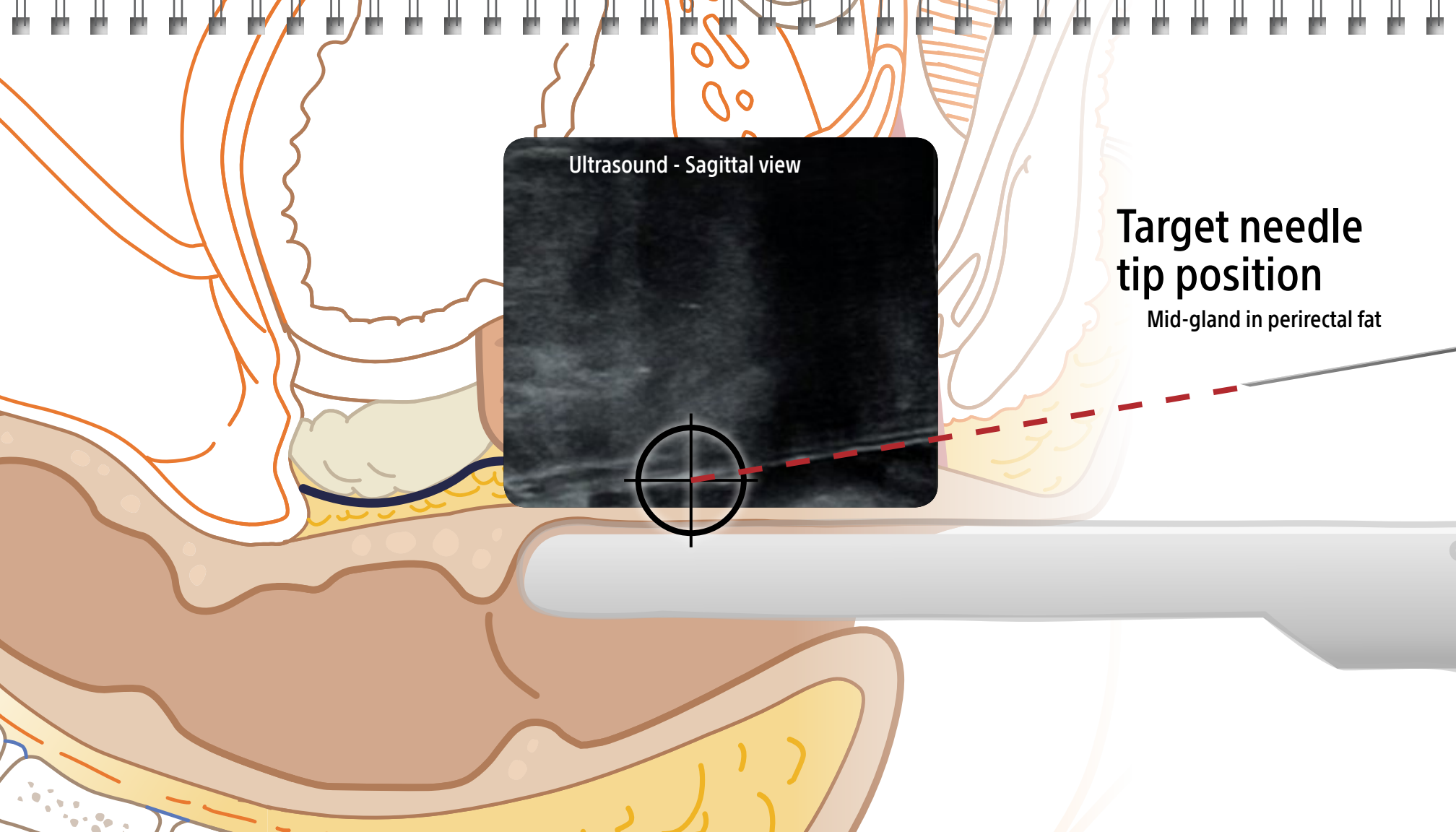


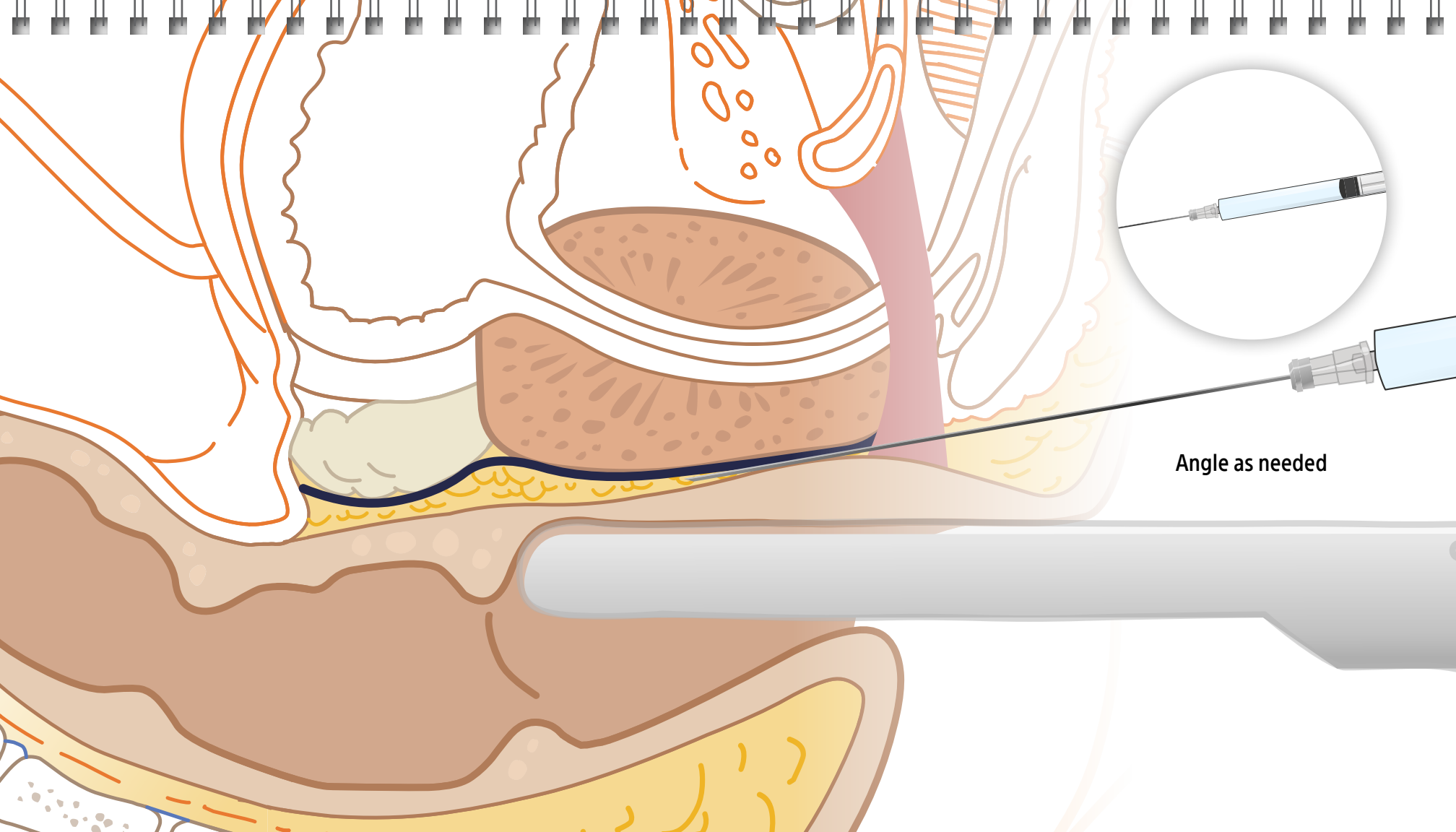


Angle as needed

Ultrasound - Axial view

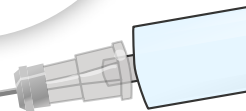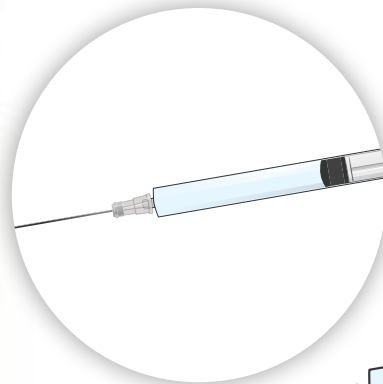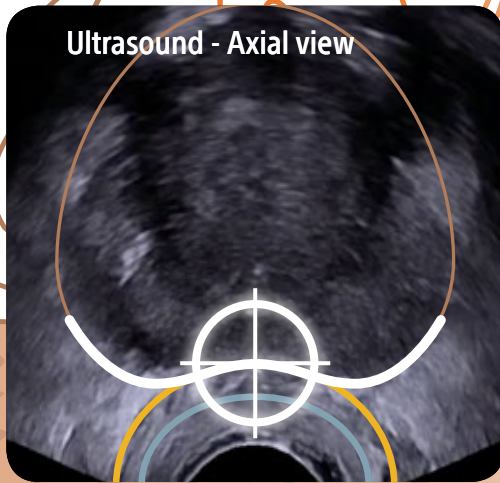


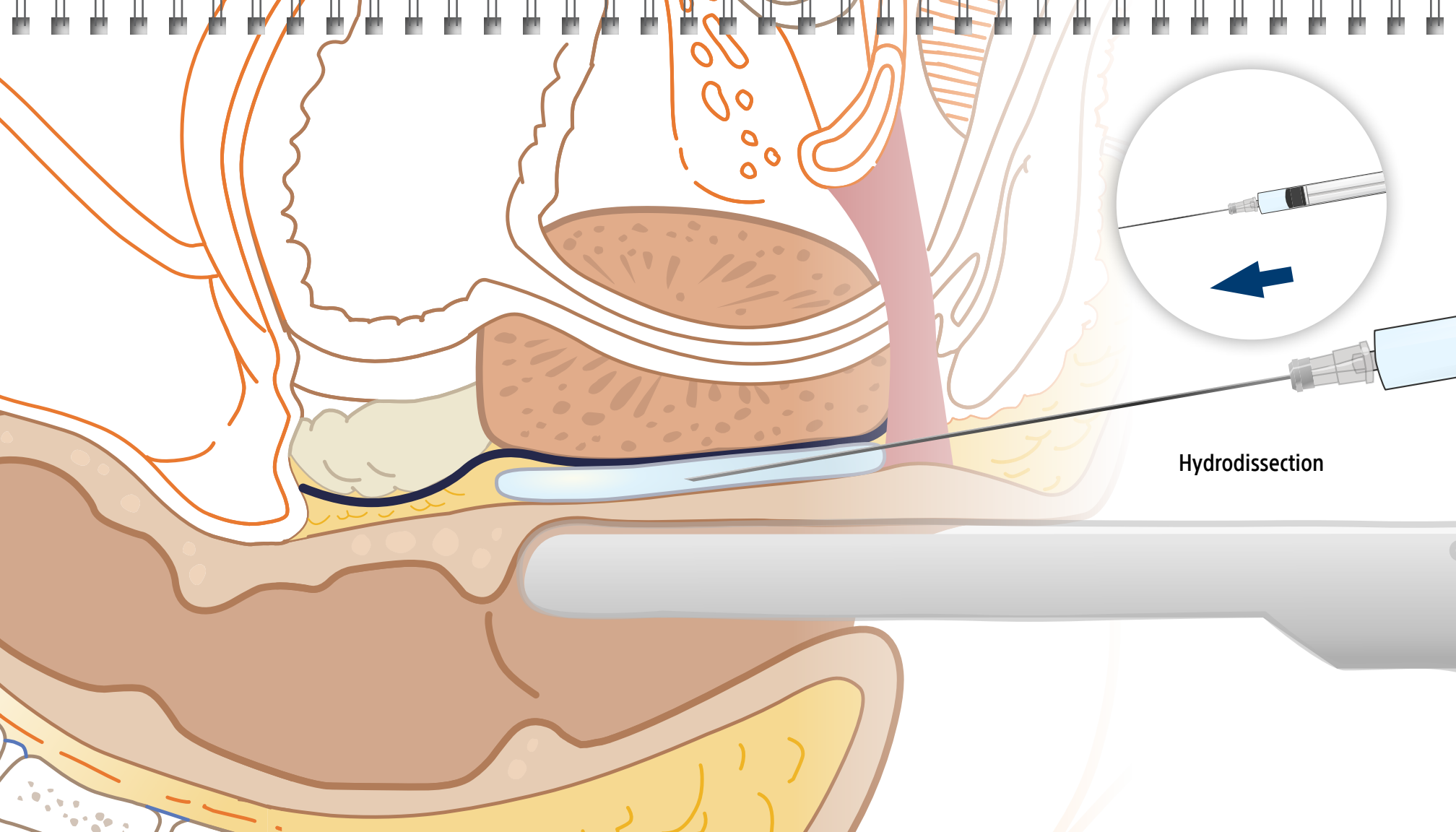


Hydrodissection

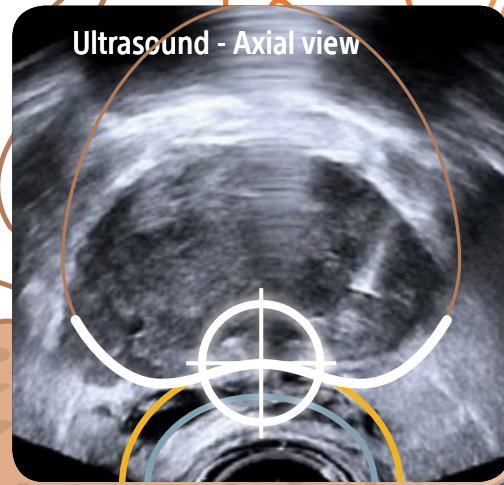


Ultrasound - Axial view

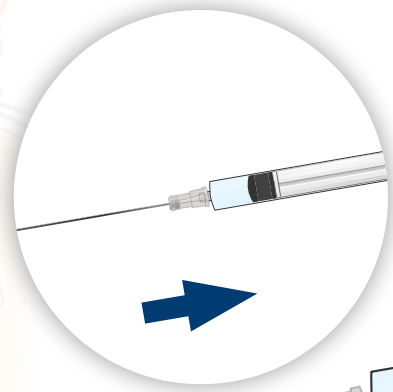


Confirm needle placement on axial view then aspirate to confirm needle is not intravascular

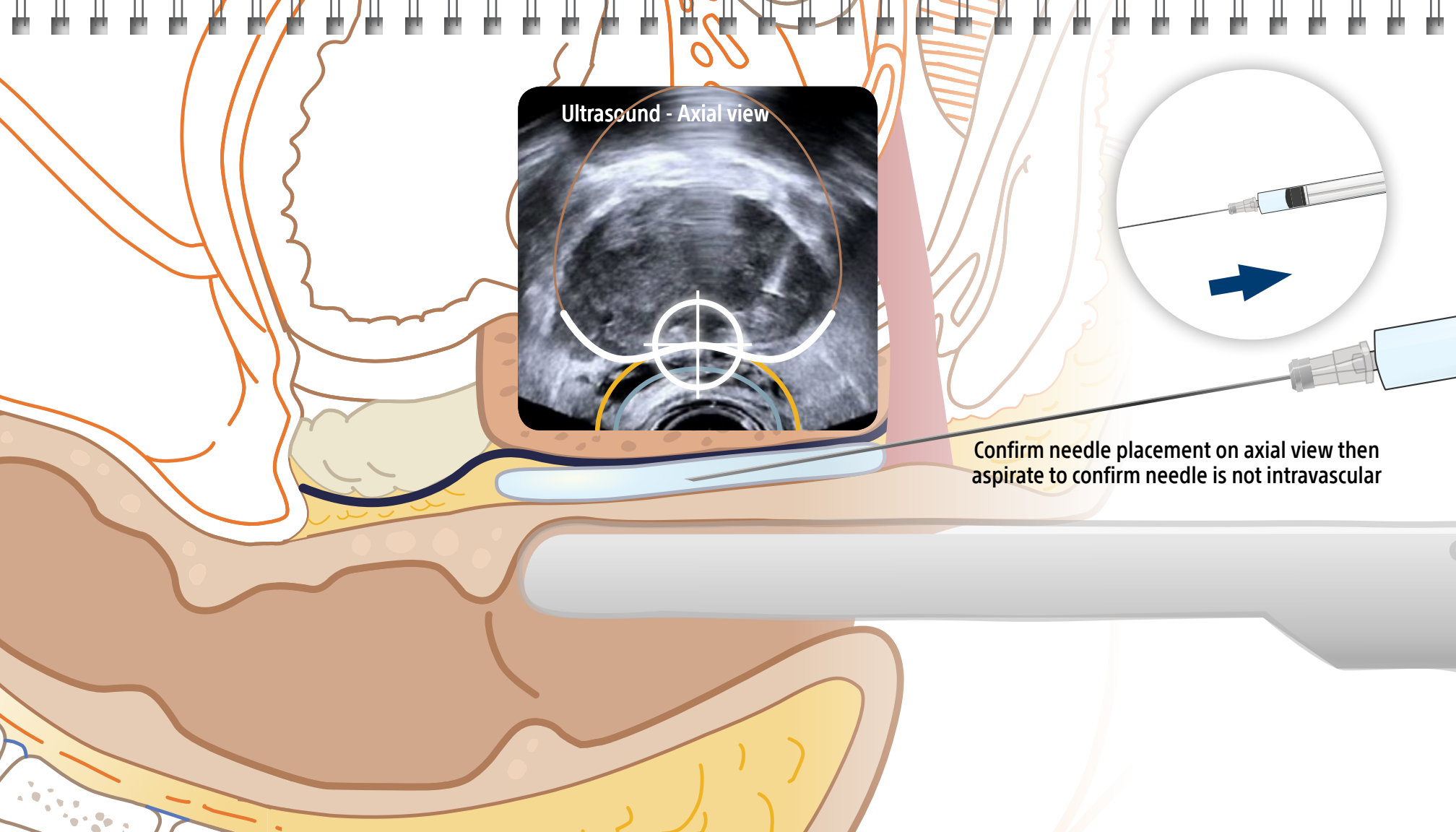

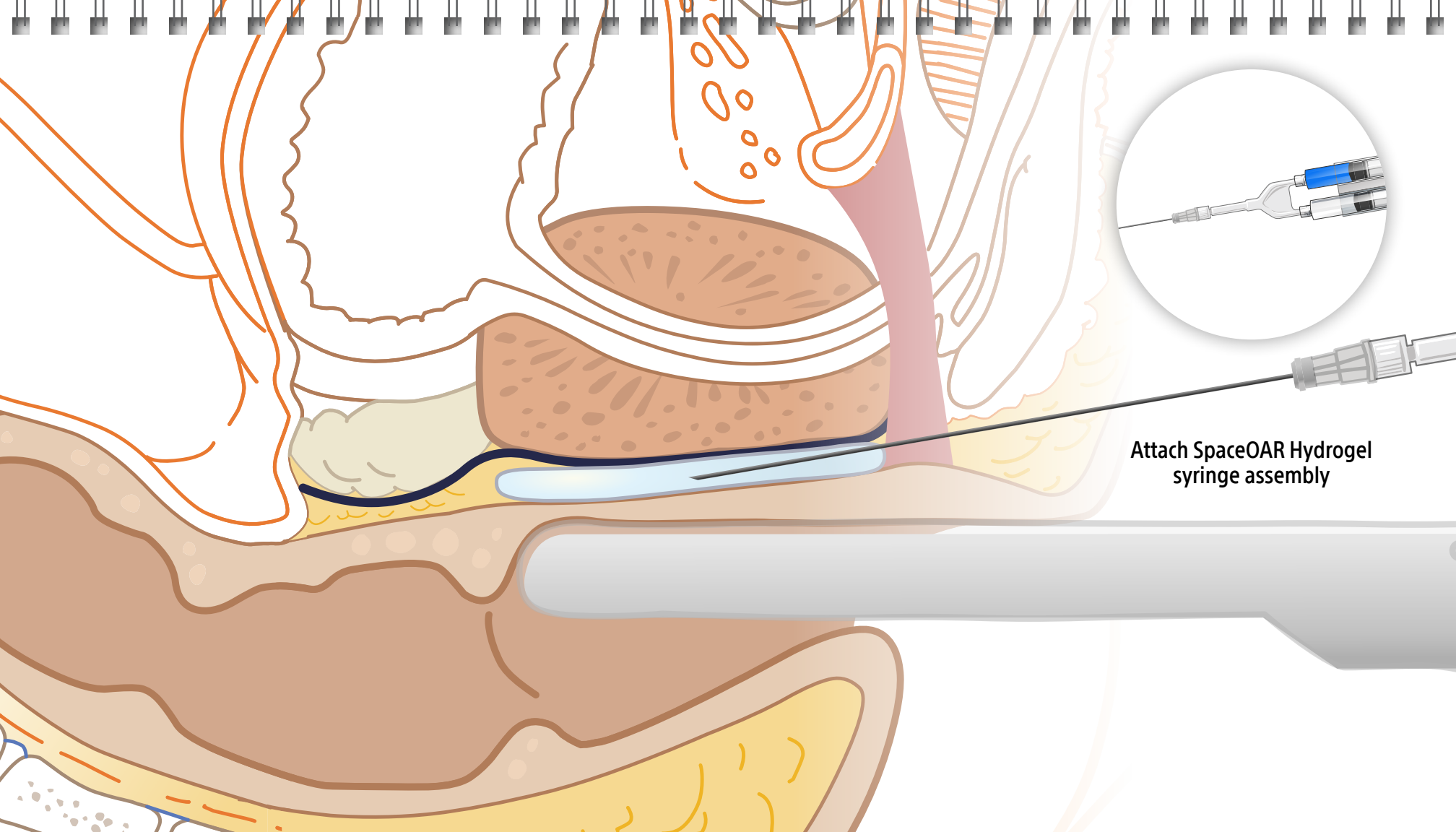


Attach SpaceOAR Hydrogel
syringe assembly

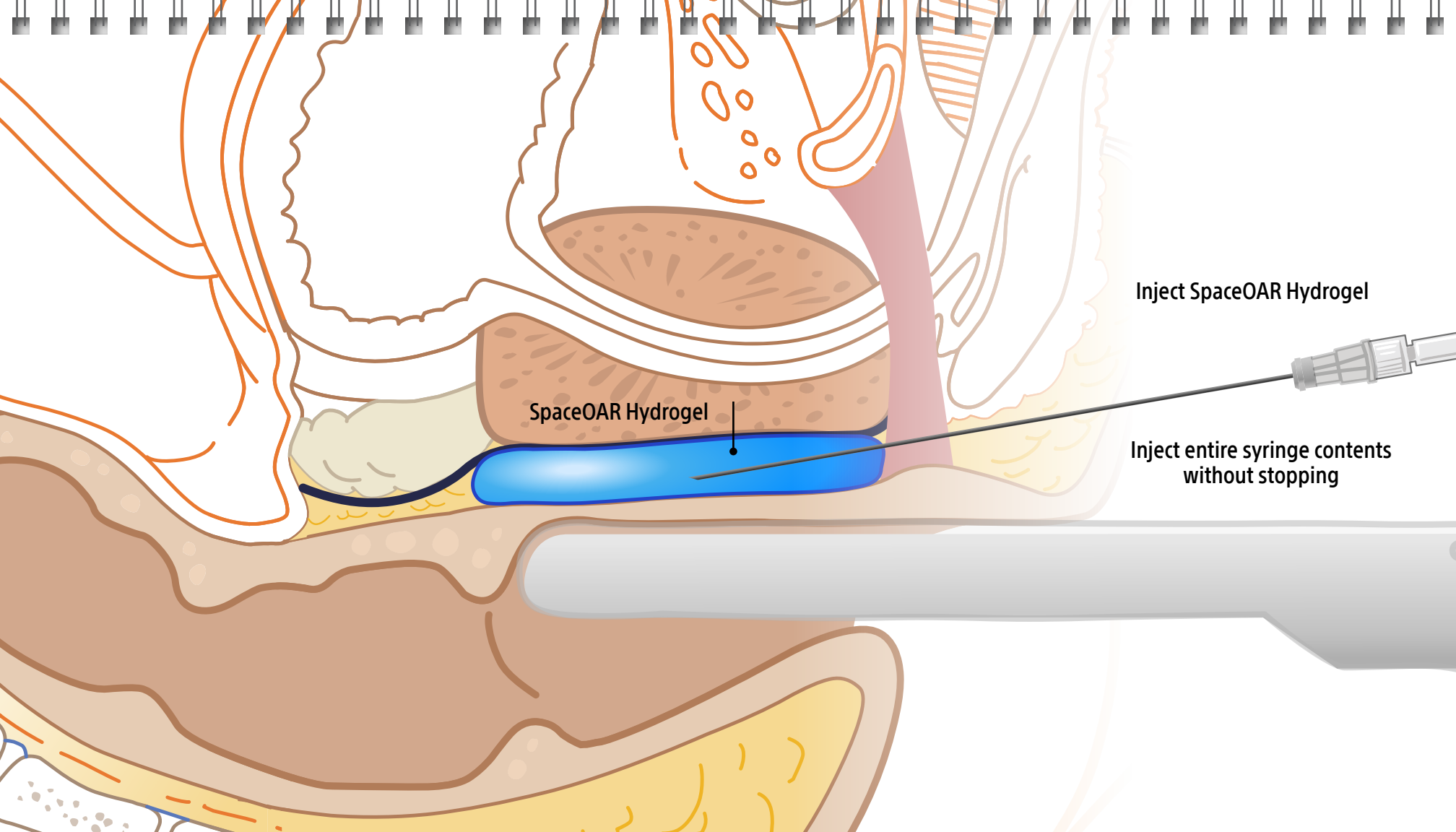


Inject SpaceOAR Hydrogel

SpaceOAR Hydrogel

Inject entire syringe contents
without stopping

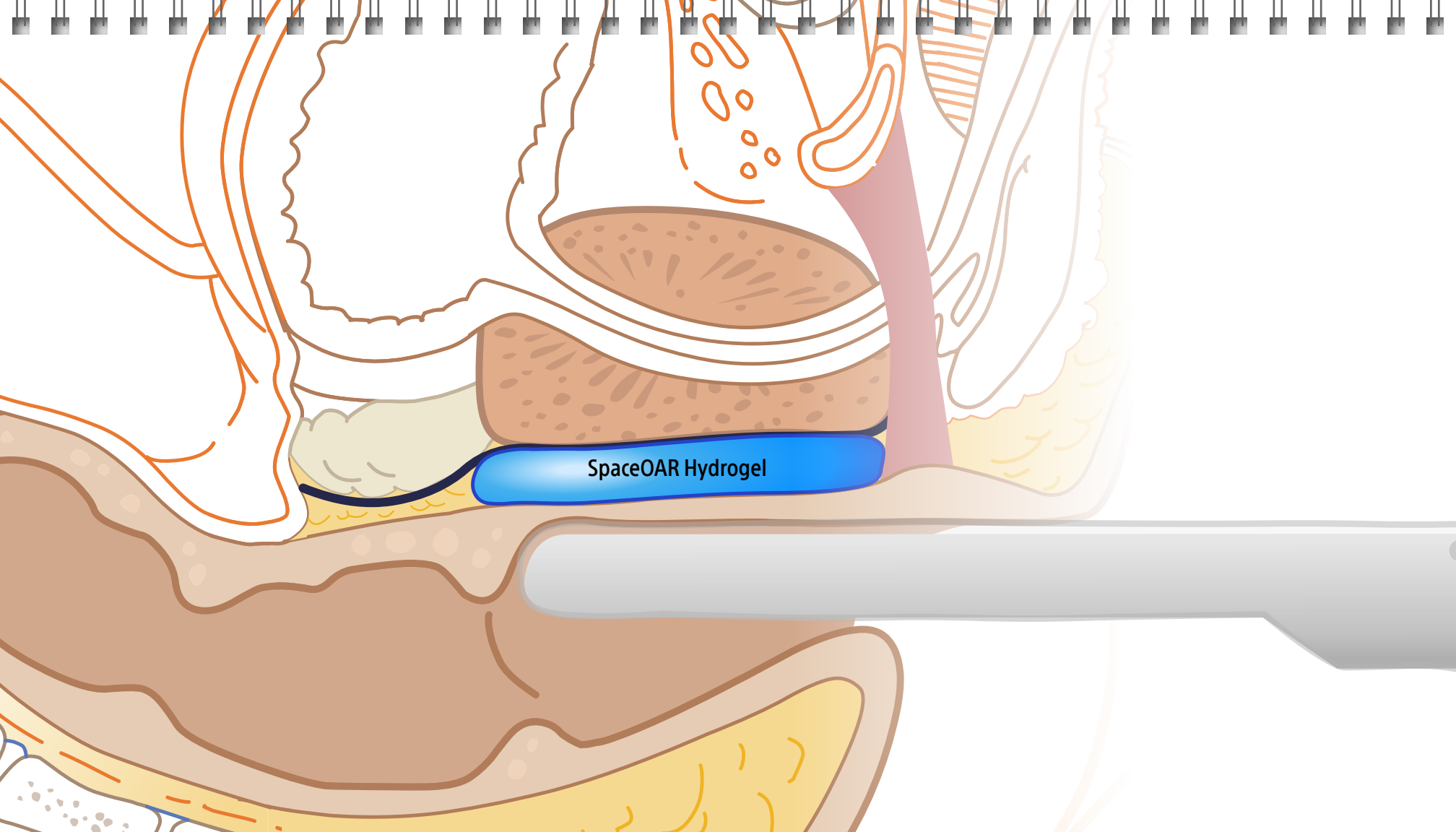

SpaceOAR Hydrogel

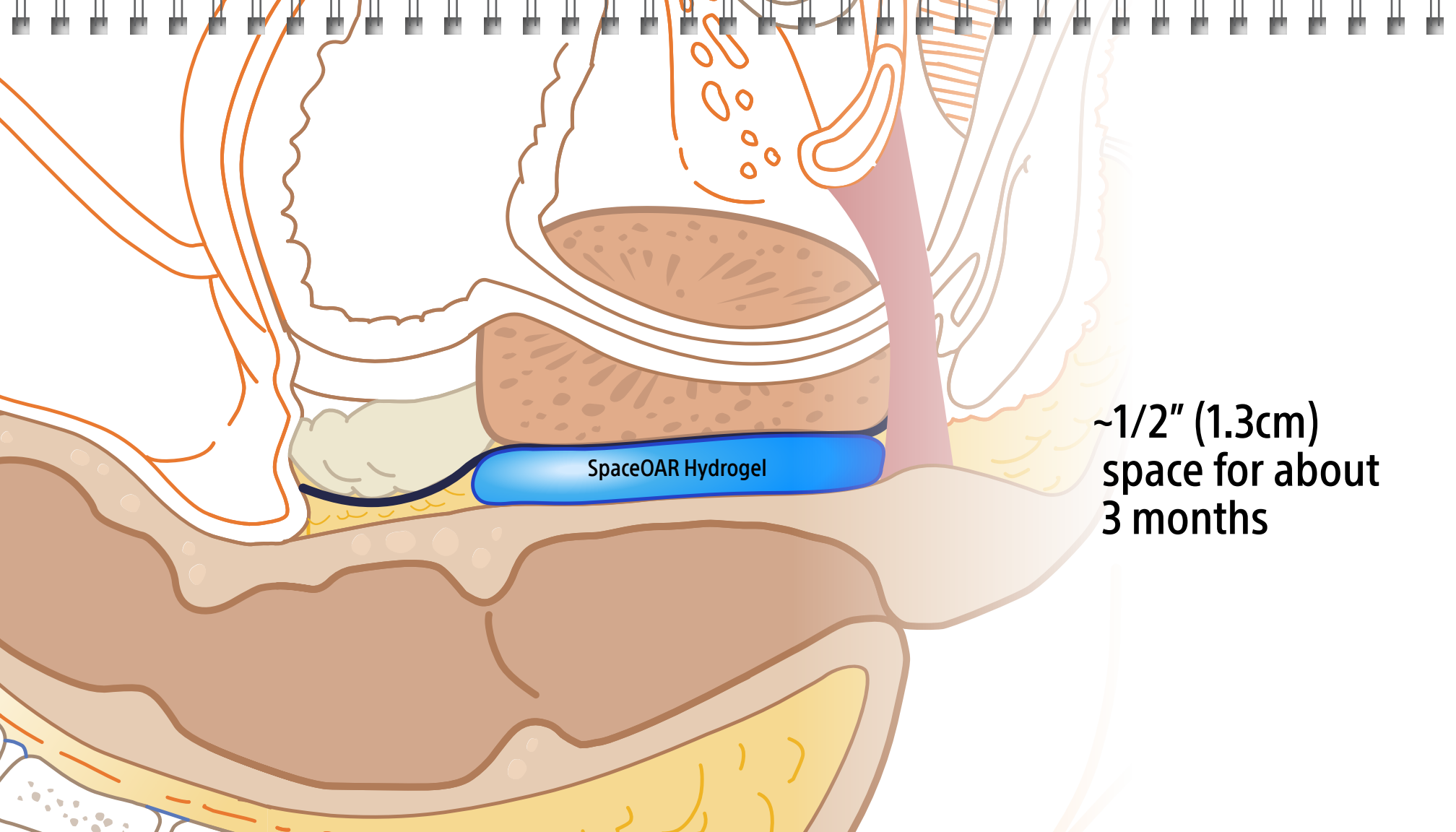

SpaceOAR Hydrogel

~1/2" (1.3cm)
space for about
3 months

As with any medical treatment, there are some risks involved with the use of SpaceOAR Hydrogel. Potential complications associated with SpaceOAR Hydrogel include, but are not limited to, pain associated with SpaceOAR Hydrogel injection; pain or discomfort associated with SpaceOAR Hydrogel, needle penetration of the bladder, prostate, rectal wall, rectum, or urethra; injection of SpaceOAR Hydrogel into the bladder, prostate, rectal wall, rectum, or urethra; local inflammatory reactions; infection; injection of air, fluid or SpaceOAR Hydrogel intravascularly; urinary retention; rectal mucosal damage, ulcers, necrosis; bleeding; constipation; and rectal urgency.

Caution: U.S. Federal law restricts this device to sale by or on the order of a physician.

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for use only in countries with applicable health authority registrations. This material not intended for use in France.

Products shown for INFORMATION purposes only and may not be approved or for sale in certain countries. Please check availability with your local sales representative or customer service.

All images are the property of Boston Scientific. All trademarks are the property of their respective owners.

**Boston
Scientific**
Advancing science for life™

Boston Scientific Corporation
300 Boston Scientific Way
Marlborough, MA 01752
www.bostonscientific.com

© 2019 Boston Scientific Corporation
or its affiliates. All rights reserved.

UROPH-620107-AB NOV 2019